



# Coronavirus (2019-nCoV)-Antigentest - Summary Data

The image shows the front of a product box for the Hotgen Coronavirus (2019-nCoV)-Antigentest. The box is white with a green and yellow gradient at the bottom. It features the Hotgen logo at the top left, a butterfly icon, and a circular illustration of a woman using the nasal swab test. The text is multilingual, providing instructions and component lists in English, German, and Spanish. A QR code for an operation video is located at the bottom right. Regulatory symbols like IVD, CE, and REF are also present.

**Hotgen**

**For self-testing use**

**Coronavirus (2019-nCoV)-Antigentest-**

Coronavirus (2019-nCoV)-Antigentest- Pour une utilisation d'auto-test  
Coronavirus (2019-nCoV)-Antigentest- Für den Selbsttest  
Coronavirus (2019-nCoV)- Test antigenico rapido- Per uso autodiagnostico  
Coronavirus (2019-nCoV)- Prueba de Antígeno- Para uso de autodiagnóstico

**Components**  
SWS-Cat-2 antigen test cassette 1 test  
Sample and action buffer 0.3 mL x 1 pcs  
Disposable Swab 1 test  
Biohazard specimen bag 1 pcs  
Instructions for use 1 pcs

**Componentes**  
Cassette de test d'antigène SARS-CoV-2 1 test  
Sérum d'action et tampon 0.3 mL x 1 PCS  
Écouvillon stérile Swab 1 test  
Sac d'échantillon Biohazard 1 PCS  
Mode d'emploi 1 PCS

**Bestandteile**  
SWS-Cat-2 Antigen-Testkassette 1 Test  
Nährchen mit Extraktionslösung 1 Stk  
Steriler Einwegswab 1 Test  
Müllbeutel für Biohazard 1 Stk  
Gebrauchsanweisung 1 Stk

**Componenti**  
Cart per test antigeno SARS-CoV-2 1 test  
Provetta di estrazione con buffer 1 pcs  
Tempone sterile monouso 1 test  
Sacchetto per rifiuti a rischio biologico 1 pcs  
Istruzioni per l'uso 1 pcs

**Componentes**  
Cassete de proba de antigén SARS-CoV-2 1 prueba  
Tampón de extracción de muestra 1 pcs  
Hisopo estéril desechable 1 test  
Bolsa para muestras de riesgo biológico 1 pcs  
Instrucciones de uso 1 pcs

IVD

6 970287 534073 **CE** 0123 **REF** HCCG134S0101

Operation Video

Beijing Hotgen Biotech Co., Ltd.

# 目录

(Contents)

1、产品彩页 (Product Brochure) .....	1-2
2、公司介绍 (Company Profile) .....	3
3、产品照片 (Product Photos) .....	4
4、包装信息 (Packing Information) . .....	5
5、热景CE自测证 (HOTGEN CE self test certificate) .....	6
6、ISO13485 认证 (ISO13485 Certificate) .....	8

# Company Profile


Beijing Hotgen Biotech Co., Ltd. (abbreviated as Hotgen Biotech, stock code: 688068) was established in June 2005, which is a high-tech enterprise focusing on the research& development, manufacture and sales of medical and public safety inspection products of in vitro diagnostics (IVD) in the field of biomedicine, as well as landed on the China Sci-Tech innovation board (STAR Market) in September 2019.

After several years of Research& development, Hotgen Biotech has developed an in vitro diagnostic reagent bioactive raw material development platform, a sugar chain abnormal protein detection (sugar capture) R&D technology platform, a Magnetic particles chemiluminescence R&D technology platform, a Up-converting Phosphor R&D technology platform, and a colloidal gold immune layer, The eight major technology platforms, such as the precipitation R&D platform, enzyme-linked immunoassay R&D technology platform, molecular diagnostics R&D platform, and instrument R&D technology platform, form a closed-loop system for in vitro diagnostic R&D and production. Hotgen Biotech has established a complete full level immunodiagnostic technology platform, from high-precision Up-converting Phosphor POCT (UPT series) to small, medium and large single- cartridge chemiluminescence platforms (MQ60 series), and then to large-scale full-automatic chemiluminescence Platform (C2000), which realizes the application of the immune diagnostic platform in the field of full diagnostic scenarios. Supporting products are widely used in the clinical and public safety fields. Specific users include hospitals at all levels, township health centers, third-party testing centers, and medical institutions, as well as medical and health institutions, as well as disease control centers, public security, fire protection, military, ports, food and medicine. Supervision, food and feed enterprises and other public safety fields.

Hotgen Biotech has won the second prize of the National Technology Invention Award, the Gold Medal of Independent Innovation, and the second prize of the Chinese Medical Science and Technology Award; In 2018, Hotgen Biotech was awarded the second prize of the "Technical Invention Category of China Rare Earth Science and Technology Award" by the China Rare Earth Society; Top 100 Private Scientific and Technological Innovations "and" Top 100 Medical Enterprises of the Future "; and" Postdoctoral Scientific Research Workstation "; major science and technology projects in the 12th and 13th five years, 863 plan, science and technology projects of the Beijing Science and Technology Commission, and Zhongguancun High Precision The project's major cutting-edge original technological achievements transformation and industrialization projects.

In the face of the COVID-19 epidemic situation, Beijing Hotgen Biotech Co.,Ltd has organized R&D developed a variety of Covid-19 detection reagents, including Novel Coronavirus 2019-nCoV Antigen Test (Colloidal Gold), Novel Coronavirus 2019-nCoV Antigen Test (Up-converting Phosphor Immunochromatographic Technology), Coronavirus disease(COVID-19) Antibody Test (Colloidal Gold), Coronavirus disease(COVID-19) IgM/IgG Antibody Rapid Test (Colloidal Gold), Novel Coronavirus 2019-nCoV Antibody Test (Up-converting Phosphor Immunochromatographic Technology), Novel Coronavirus (SARS-CoV-2) Neutralizing Antibodies Test (Colloidal Gold), Coronavirus disease( COVID-19) Nucleic Acid Test Kit (PCR-Fluorescent Probe Method), Disposable virus sampling tube, Nucleic acid Automatic Purification System, Nucleic acid extraction reagent, Biological Sample Releaser kit, etc.It is imperative to fight the epidemic Helping the global fight against epidemics!

Since its establishment, the company has continuously grown its business and has now achieved group development. At present, Hotgen (Langfang), Hotgen (Jilin), Weikekang Technology, Shunjing Biological and many other subsidiaries have been established. Hotgen Biotech marketing and service network has covered all regions of the country. Each province is equipped with professional technical service engineers, academic engineers, etc. who are responsible for pre-sales and after-sales work to meet customer needs. The company takes "developing biotechnology and benefiting human health" as its mission, "quality determines the company's life and death, customers determine the company's success or failure, talents determine the company's rise and fall, innovation determines the company's future" as its core values, and "tests because of me advanced" as its philosophy , High ambitions, technological entrepreneurship, and industrial prosperity!

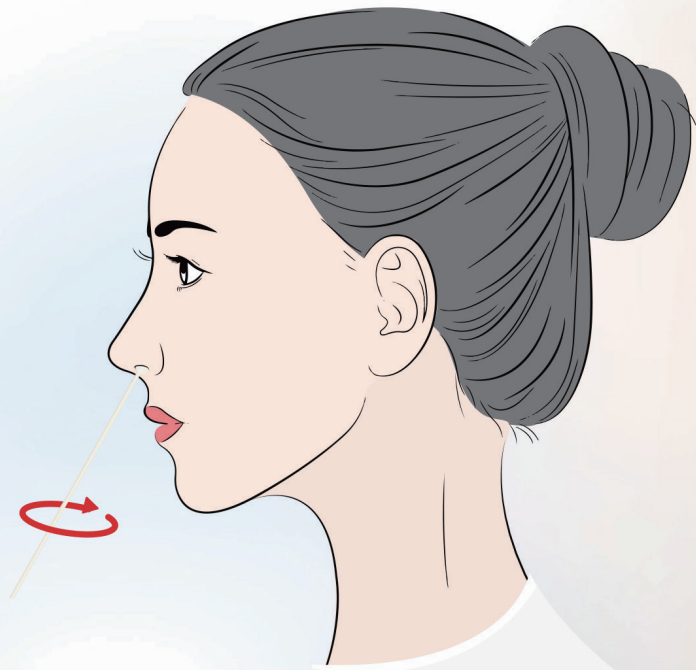

## Coronavirus (2019-nCoV)-Antigentest-



### Product Features

- High Accuracy, Specificity and Sensitivity ●
- No need instrument, get results in 15 minutes ●
- Room temperature storage ●
- Sample : Human Anterior Nares Swab ●
- Detect the presence of viral proteins ●
- Identify acute or early infection ●

# Specimen Requirements and Test Procedure

## Specimen Requirements

1

Please wash and dry your hands thoroughly before the test.




2

Read the Instructions for Use carefully.

3



Please scan the QR code to watch the operation video.

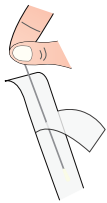

4

Take the swab out of the package and do not touch the sampling end.


5



Carefully insert the swab 1.5 cm into the nostril until resistance is felt. Under moderate pressure, swab the surface of the nostril 4 to 6 times for at least 15 seconds.

6

Repeat sampling with the same swab in the other nostril.

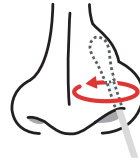


7

Open the larger end of the sampling tube, and insert the swab after collecting the sample into the sampling tube.

8

Soak the swab in the buffer for at least 15 seconds, stir the swab several times, and squeeze the swab head 3 times.

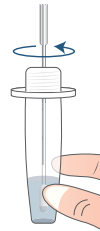

9

When remove the swab, squeeze the sides of the sampling tube.



10

Close the sampling tube with the tube lid




11

Open the foil pouch and place the test cassette on a flat surface.



12

Open the tube lid at the front end of the sampling tube, and drop 4 drops of the treated sample into the sample well (S) of the test cassette.

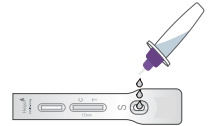

13



After reacting at room temperature for 15 minutes, observe the results. A result after 30 minutes is invalid.

14

After the test, pack all the components of this test in the Biohazard specimen bag for contaminated waste and dispose of this bag closed with the residual waste. Not reusable.




15




Wash or disinfect hands again.


## Interpretation of result

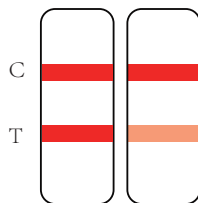

Positive

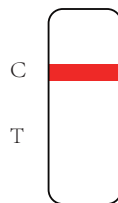

Negative

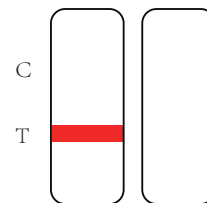

Invalid

## Product information

Product name	Test samples	Specification	Storage conditions
Coronavirus (2019-nCoV)-Antigentest-	Human Anterior Nares Swab	1T/kit	4-30 C

# Coronavirus (2019-nCoV)-Antigentest - Product Photos

**Hotgen**

**For self-testing use**

**Coronavirus (2019-nCoV)-Antigentest-**

Coronavirus (2019-nCoV)-Antigentest- Pour une utilisation d'auto-test  
 Coronavirus (2019-nCoV)-Antigentest- Für den Selbsttest  
 Coronavirus (2019-nCoV)- Test antigenico rapido- Per uso autodiagnostico  
 Coronavirus (2019-nCoV)- Prueba de Antígeno- Para uso de autodiagnóstico

**Components**  
 SARS-CoV-2 antigen test cassette 1 test  
 Sample extraction buffer 0.3 mL x 1 pcs  
 Disposable Sterile Swab 1 piece  
 Release/operation bag 1 pcs  
 Instructions for use 1 pcs

**Composantes**  
 Cassette de test d'antigène SARS-CoV-2 1 test  
 Tempois d'extraction d'échantillon 0.3 mL x 1 PCS  
 Écouvillon stérile jetable 1 pièce  
 Sac d'isolation/Release/operation 1 PCS  
 Mode d'emploi 1 PCS

**Bestandteile**  
 SARS-CoV-2 Antigen Testkassette 1 Test  
 Röhrechen mit Extraktionslösung 1 Stk  
 Steriler Einwegspatel 1 Test  
 Müllbeutel für Geländeabfall 1 Stk  
 Gebrauchsanweisung 1 Stk

**Componenti**  
 Cass per test antigenico SARS-CoV-2 1 test  
 Provanta di estrazione con buffer 1 pcs  
 Tampone sterile monouso 1 test  
 Sacchetto per rifiuti a rischio biologico 1 pcs  
 Istruzioni per l'uso 1 pcs

**Componentes**  
 Cassete de pruebas de antígeno SARS-CoV-2 1 prueba  
 Tiempo de extracción de muestras 1 pcs  
 Hisopo estéril desechable 1 prueba  
 Bolsa para muestras de riesgo biológico 1 pcs  
 Instrucciones de uso 1 pcs

IVDI, CE 0123, REF: HGCG134SD101, 6 97 0297 534 073, Operation video

**Test Procedure**

**1 Sample collection:**

- Gently insert the entire soft tip of the swab into one nostril for 1.5cm until you feel a bit of resistance.
- Using medium pressure, rub the swab slowly in a circular motion around the inside wall of your nostril 4 times for a total time of 15 seconds.
- Repeat the same process with the same swab in the other nostril.

**2** The swab after sampling is soaked below the liquid level of the sampling tube, rotated and pressing 3 times, the swab soaking time is not less than 15s, the swab head is pressed, then taken out the swab and tighten the sampling tube. The liquid in the tube is the sample after treatment.

**3** Add 4 drops of the treated sample into the sample well of the test cassette.

**4** 15 min. Observe results after 15 minutes. Result got after 30 minutes is invalid.

**5 Interpretation of result**

Positive: C (red line), T (red line)  
 Negative: C (red line), T (no red line)  
 Invalid: C (no red line), T (no red line)

MeiNet GmbH, Bismarckstrasse 10, 48163 Münster, Germany  
 Beijing Hotgen Biotech Co., Ltd., 9th Building, No. 9 Tianfu Street, Bomedical Base, Daxing District, Beijing, 102600, P.R.China

## 前鼻腔抗原胶体金检测试剂包装信息

### Coronavirus (2019-nCoV)-Antigentest - Packing Information

产品名称 Product name	规格/盒 Specifications	单位 Unit	单位包装毛重 Gross weight per unit package
Coronavirus (2019-nCoV)-Antigentest -	1T	袋/kit	0.031 kg/袋 0.031 kg / kit

前鼻腔抗原胶体金试剂盒出口包装箱 Coronavirus (2019-nCoV)-Antigentest - Export Packing Cartons							
包装箱 / 盒 Packing Carton/ box	长 length cm	宽 Width cm	高 height cm	每箱装袋 数量 Kit quantity per carton	单盒试剂 净重 Net weight of single kit	整箱净重 Net weight of the whole carton	抛重 Throwing weight
纸箱 carton	73	39	26	400袋 400T	0.031 公斤 0.031 kg	13.5公斤 13.5kg	12.3公斤 12.3kg

For reference only, shipping cost according to the actual weight of the airline department







Benannt durch/Designated by  
 Zentralstelle der Länder  
 für Gesundheitsschutz  
 bei Arzneimitteln und  
 Medizinprodukten  
 www.zlg.de  
 ZLG-BS-245.10.07



Product Service

# EC Certificate

EC Design-Examination Certificate  
 Directive 98/79/EC on In Vitro Diagnostic Medical Devices (IVDD), Annex III (6)  
 (Devices for self-testing)

**No. V9 089675 0006 Rev. 00**

**Model(s):** **Coronavirus (2019-nCoV)-Antigentest-**

**Facility(ies):** Beijing Hotgen Biotech Co.,Ltd  
 9th Building, No. 9 Tianfu Street, Biomedical Base, Daxing District,  
 102600 Beijing, PEOPLE'S REPUBLIC OF CHINA

Model Name:	REF number:
Coronavirus (2019-nCoV)-Antigentest-	HGCG134S0101
Coronavirus (2019-nCoV)-Antigentest-	HGCG134S0105
Coronavirus (2019-nCoV)-Antigentest-	HGCG134S0120
Coronavirus (2019-nCoV)-Antigentest-	HGCG134S0140

仅供外贸备案使用

